


脂肪酸甲酯混合物

(*cis/trans* Isomers) on MEGA-10 FAMES Column

Column: **MEGA-10 FAST - 0.20mm, 0.20µm, 30m**

Catalog Code: F-10-020-020-30

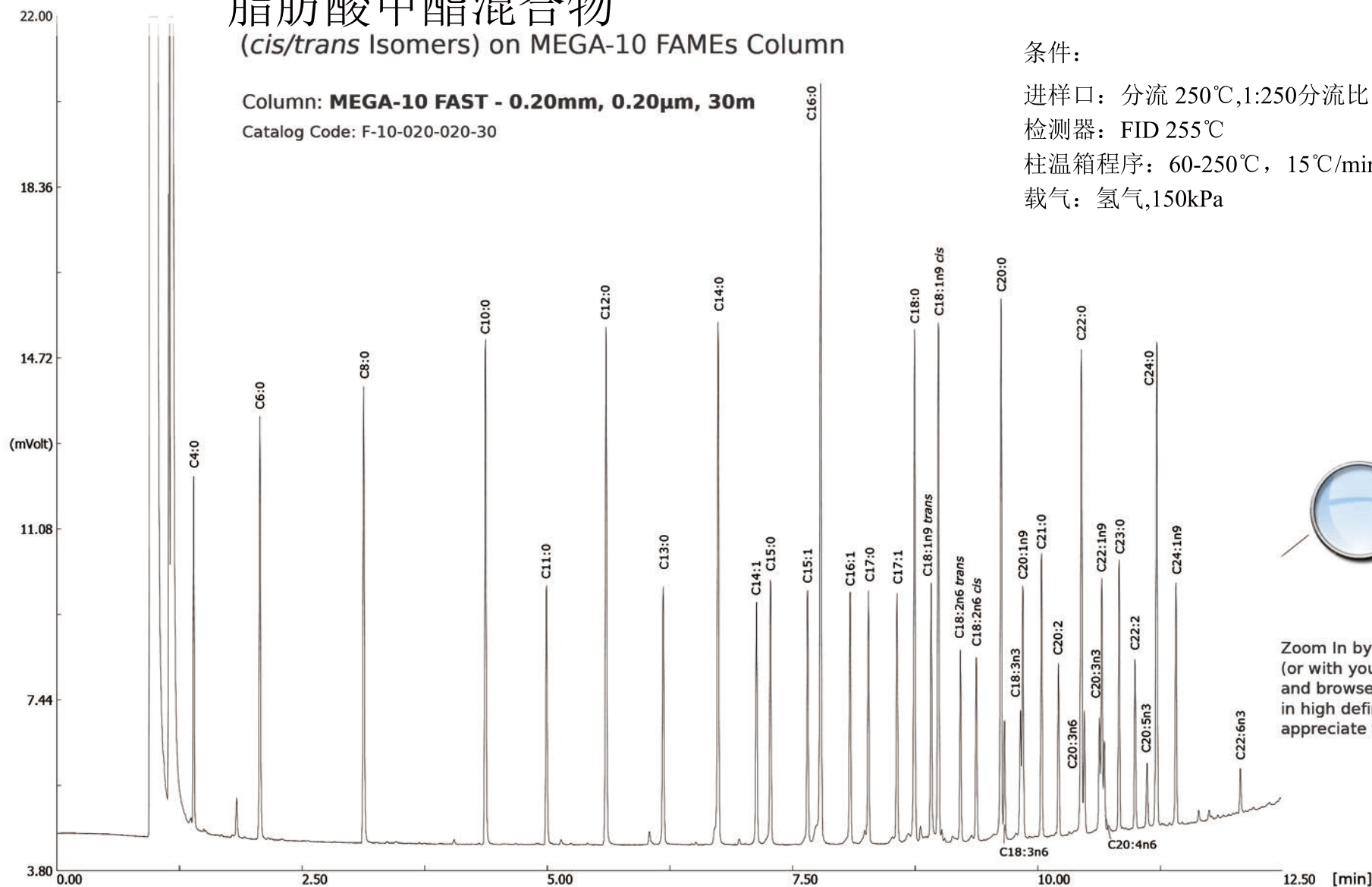

条件:

进样口: 分流 250°C, 1:250分流比, 0.5µL进样量

检测器: FID 255°C

柱温箱程序: 60-250°C, 15°C/min

载气: 氢气, 150kPa



Zoom In by clicking here
(or with your Pdf viewer zoom)
and browse the chromatogram
in high definition to better
appreciate the separation!

