


乙酰基反式环己烷(+/-)的手性分离

Column: **MEGA-DEX DMT Beta - 0.25mm, 0.25 μ m, 25m**

Catalog Code: DEX-DMT-B-025-025-25

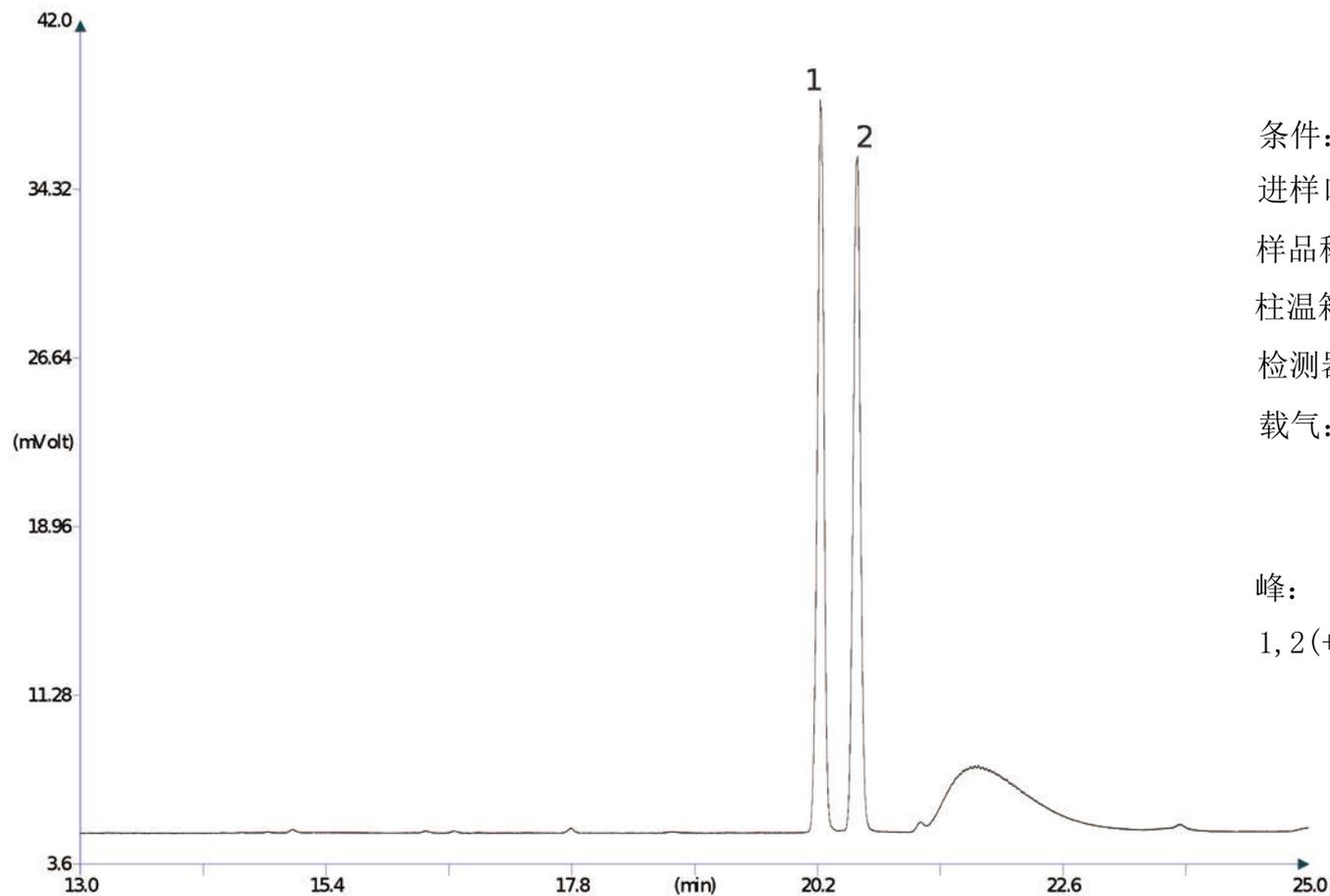

条件:

进样口: 分流250 $^{\circ}$ C, 0.5 μ L进样量, 分流比为100mL/min

样品稀释: 1mL中溶解1mg

柱温箱程序: 120 $^{\circ}$ C-200 $^{\circ}$ C, 2.5 $^{\circ}$ C/min

检测器: FID 250 $^{\circ}$ C

载气: 氦气, 110KPa

峰:

1, 2(+/-)反式环己烷