
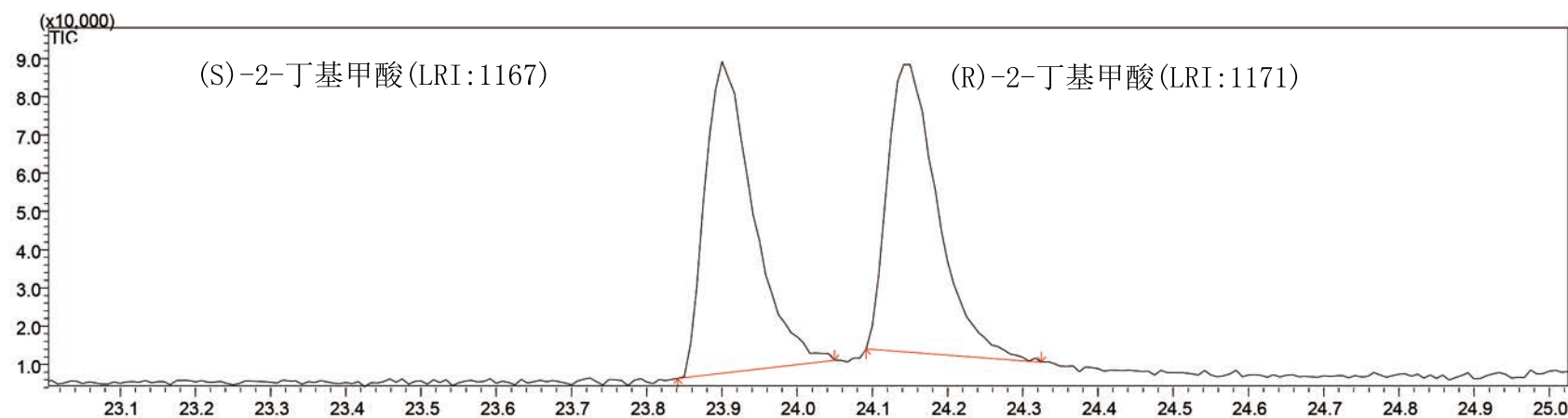


# 乙基-2-丁酸甲酯和2-甲基丁酸

GC-MS  
columns

dex xeb  
chiral columns

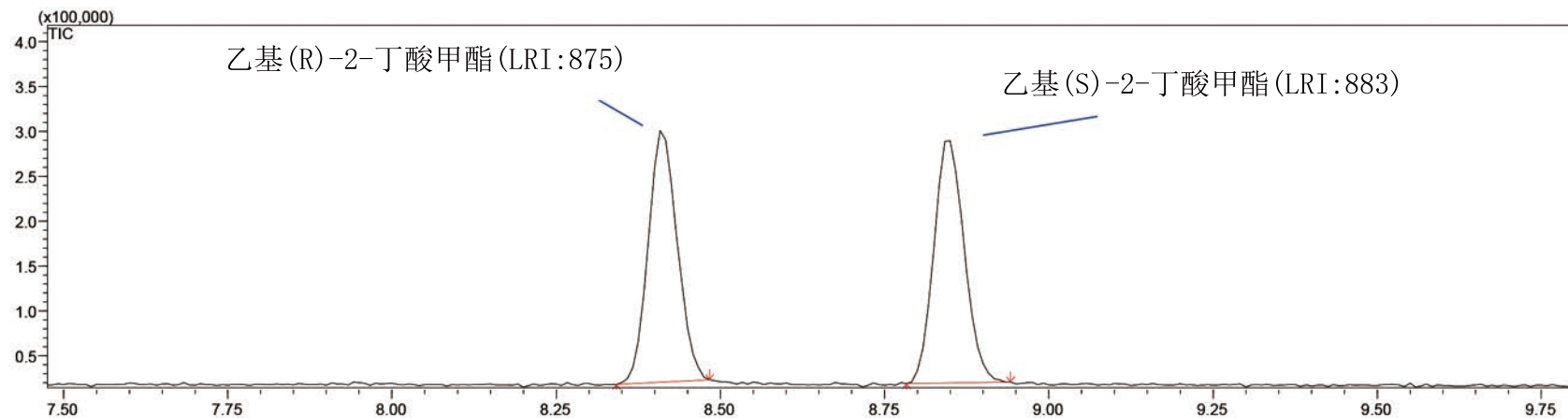

Column:

**MEGA-DEX DET Beta - 0.25mm, 0.25 $\mu$ m, 25m**

Catalog Code: DEX-DET-B-025-025-25

条件:

进样口: 分流 250 $^{\circ}$ C, 50mL/min分流比, 1 $\mu$ L进样量

柱温箱程序: 50 $^{\circ}$ C, 2 $^{\circ}$ C/min

载气: 氦气, 1mL/min

制样: 标准混合, 1mL正己烷中溶解1 $\mu$ L样品

All our MEGA-DEX chiral GC columns work perfectly with MS detectors

Acknowledgments: many thanks to Barbara Sgorbini, Cecilia Cagliero, Chiara Cordero and Carlo Bicchi, Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via P. Giuria, 9, Torino, Italy.

MEGA<sup>®</sup>

[www.mega.mi.it](http://www.mega.mi.it)