


洋甘菊精油

Column: **MEGA-WAX FAST - 0.10mm, 0.10 μ m, 5m**

Catalog Code: F-WAX-010-010-5

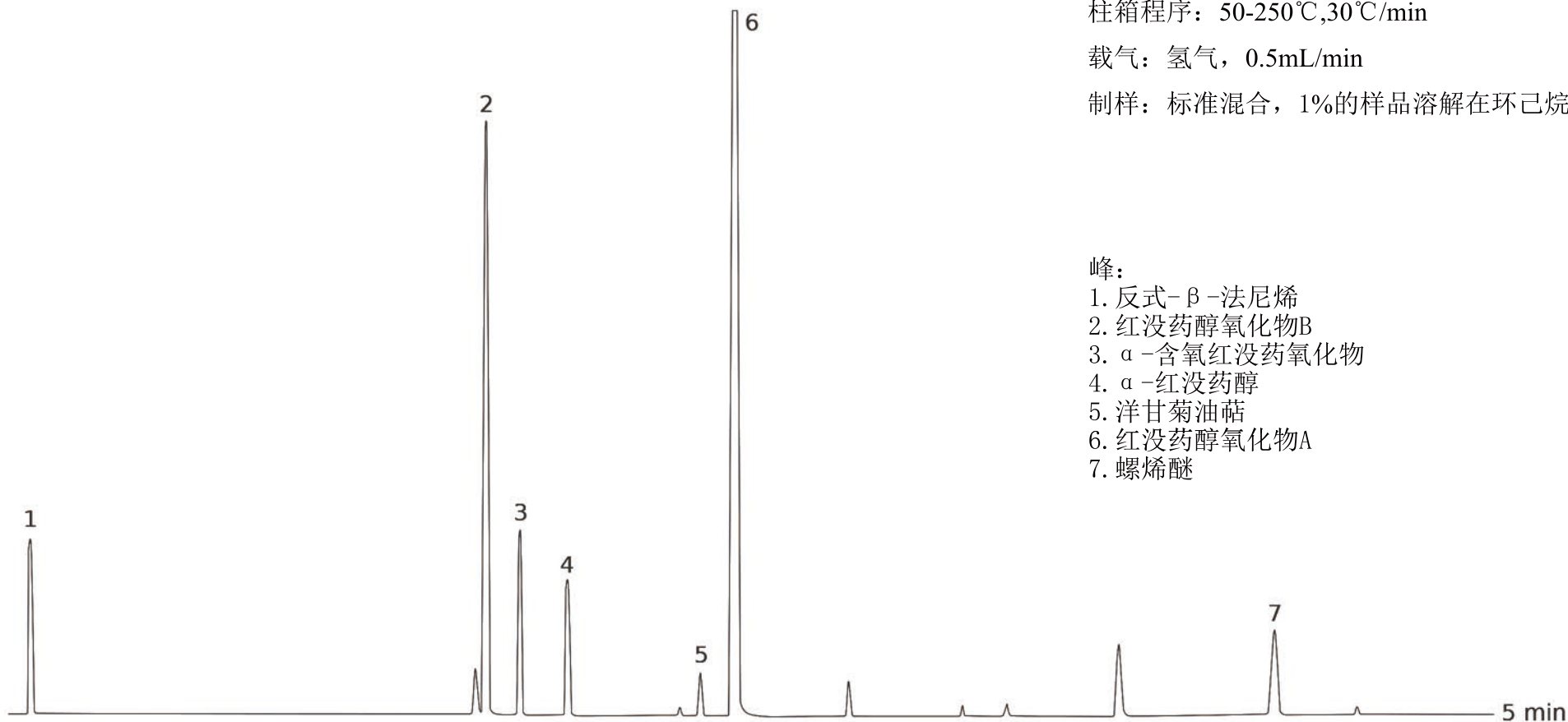

条件:

进样口: 分离230 $^{\circ}$ C,0.5 μ L进样量,1:250分流比

检测器: FID 250 $^{\circ}$ C

柱箱程序: 50-250 $^{\circ}$ C,30 $^{\circ}$ C/min

载气: 氢气, 0.5mL/min

制样: 标准混合, 1%的样品溶解在环己烷中

峰:

1. 反式- β -法尼烯
2. 红没药醇氧化物B
3. α -含氧红没药氧化物
4. α -红没药醇
5. 洋甘菊油萜
6. 红没药醇氧化物A
7. 螺烯醚

Acknowledgement: Prof. C.Bicchi, C.Brunelli et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.