


食品包装的残留溶剂

Column: **MEGA-VOC 2 FAST - 0.20mm, 1.00 μ m, 10m**

Catalog Code: F-VOC2-020-100-10

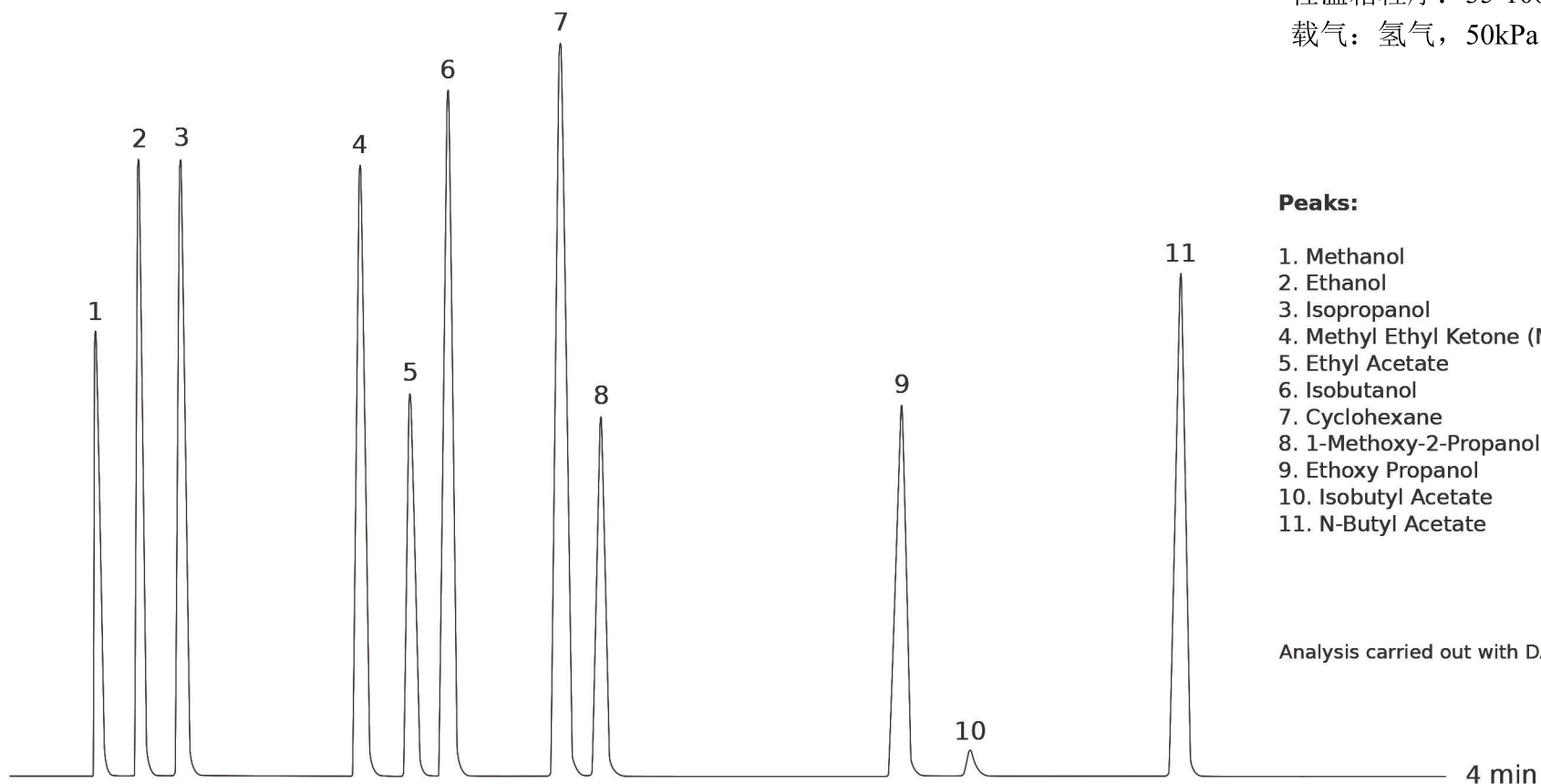

条件:

进样口: 分流 250 $^{\circ}$ C

检测器: FID 250 $^{\circ}$ C

柱温箱程序: 35-100 $^{\circ}$ C, 15 $^{\circ}$ C/min

载气: 氢气, 50kPa

Peaks:

1. Methanol
2. Ethanol
3. Isopropanol
4. Methyl Ethyl Ketone (MEK)
5. Ethyl Acetate
6. Isobutanol
7. Cyclohexane
8. 1-Methoxy-2-Propanol
9. Ethoxy Propanol
10. Isobutyl Acetate
11. N-Butyl Acetate

Analysis carried out with DANI MASTER-GC.