

# 杀虫剂

Column: **MEGA-SE54 FAST - 0.10mm, 0.10µm, 5m**

Catalog Code: F-SE54-010-010-5

条件:


进样口: 分流 230°C, 1:200分流比, 0.5µL进样量

检测器: FID 250°C

柱温箱程序: 50-250°C, 40°C/min

载气: 氢气, 0.5mL/min

制样: 标准混合, 1%正己烷中溶解

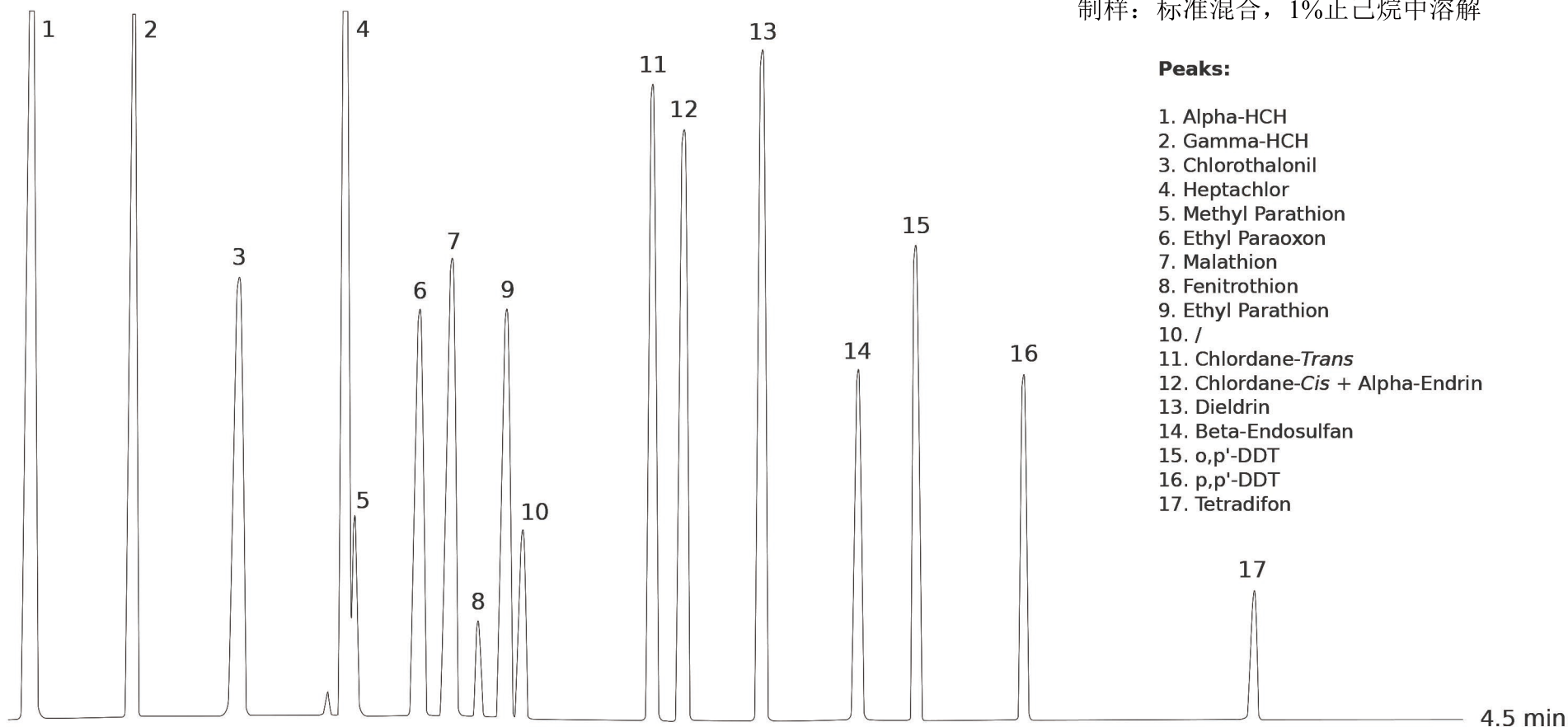

### Peaks:

1. Alpha-HCH
2. Gamma-HCH
3. Chlorothalonil
4. Heptachlor
5. Methyl Parathion
6. Ethyl Paraoxon
7. Malathion
8. Fenitrothion
9. Ethyl Parathion
10. /
11. Chlordane-Trans
12. Chlordane-Cis + Alpha-Endrin
13. Dieldrin
14. Beta-Endosulfan
15. o,p'-DDT
16. p,p'-DDT
17. Tetradifon

Acknowledgement: Prof. C.Bicchi, C.Brunelli et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.